sanitas

SBM 18

DE	Blutdruckmessgerät	
	<i>Gebrauchsanweisung</i>	2
EN	Blood pressure monitor	
	<i>Instructions for use</i>	14
FR	Tensiomètre	
	<i>Mode d'emploi</i>	25
IT	Misuratore di pressione	
	<i>Istruzioni per l'uso</i>	37
RU	Прибор для измерения кровяного давления в плечевой артерии	
	<i>Инструкция по применению</i>	48
PL	Ciśnieniomierz na przedramię	
	<i>Instrukcja obsługi</i>	61
NL	Bloeddrukmeter voor de bovenarm	
	<i>Gebruiksaanwijzing</i>	73

Inhalt

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Contents

1. Getting to know your device	14
2. Important notes	14
3. Device description	18
4. Preparing the measurement	19
5. Measuring blood pressure	19
6. Evaluating results	21
7. Saving, displaying and deleting measured values	22
8. Error messages/troubleshooting	23
9. Cleaning and storing the device and cuff	23
10. Technical specifications	23
11. Replacement parts and wearing parts	24
12. Warranty / Service	24

Included in delivery

- Blood pressure monitor
- Upper arm cuff
- 4 x 1.5 V LR6 AA batteries
- Storage bag
- Instructions for use

Please read these instructions for use carefully and keep them for later use, be sure to make them accessible to other users and observe the information they contain.

1. Getting to know your device

Check that the device packaging has not been tampered with and make sure that all contents are present. Before use, ensure that there is no visible damage to the device or accessories and that all packaging material has been removed. If you have any doubts, do not use the device and contact your retailer or the specified Customer Services address.

The upper arm blood pressure monitor is used to carry out non-invasive measurement and monitoring of the arterial blood pressure values in adults.

This allows you to quickly and easily measure your blood pressure, save the measured values and display the development and average values of the measured values taken.

You are also warned of possible existing cardiac arrhythmia.

The recorded values are classified and evaluated graphically. Store these instructions for use for future reference and make them accessible to other users.

2. Important notes

Signs and symbols

The following symbols are used in these instructions for use, on the packaging and on the type plate for the device and the accessories:

Attention

	Note Note on important information
	Observe the instructions for use
	Application part, type BF
	Direct current
	Disposal in accordance with the Waste Electrical and Electronic Equipment EC Directive – WEEE
	Dispose of packaging in an environmentally friendly manner
	Manufacturer
	Permissible storage and transport temperature and humidity
	Permissible operating temperature and humidity
IP21	Protected against solid foreign objects 12.5 mm in diameter and larger, and against vertically falling drops of water
	Serial number

 0483	CE labelling This product satisfies the requirements of the applicable European and national directives.
---	--

Notes on use

- In order to ensure comparable values, always measure your blood pressure at the same time of day.
- Do not take a measurement within 30 minutes of eating, drinking, smoking or exercising.
- Before the initial blood pressure measurement, make sure always to rest for about 5 minutes.
- Furthermore, if you want to take several measurements in succession, make sure always to wait for at least 1 minute between the individual measurements.
- Repeat the measurement if you are unsure of the measured value.
- The measured values taken by you are for your information only – they are no substitute for a medical examination. Discuss the measured values with your doctor and never base any medical decisions on them (e.g. medicines and their administration).
- Using the blood pressure monitor outside your home environment or whilst on the move (e.g. whilst travelling in a car, ambulance or helicopter, or whilst undertaking physical activity such as playing sport) can influence the measurement accuracy and cause incorrect measurements.

- This device is not intended for use by people (including children) with restricted physical, sensory or mental skills or a lack of experience and/or a lack of knowledge, unless they are supervised by a person who has responsibility for their safety or they receive instructions from this person on how to use the device. Supervise children around the device to ensure they do not play with it.
- Do not use the blood pressure monitor on newborns or patients with pre-eclampsia. We recommend consulting a doctor before using the blood pressure monitor during pregnancy.
- Cardiovascular diseases may lead to incorrect measurements or have a detrimental effect on measurement accuracy. The same also applies to very low blood pressure, diabetes, circulatory disorders and arrhythmias as well as chills or shaking.
- The blood pressure monitor must not be used in connection with a high-frequency surgical unit.
- Only use the device on people who have the specified upper arm measurement for the device.
- Please note that when inflating, the functions of the limb in question may be impaired.
- During the blood pressure measurement, the blood circulation must not be stopped for an unnecessarily long time. If the device malfunctions remove the cuff from the arm.
- Avoid any mechanical restriction, compression or bending of the cuff line.
- Do not allow sustained pressure in the cuff or frequent measurements. The resulting restriction of the blood flow may cause injury.
- Make sure that the cuff is not placed on an arm in which the arteries or veins are undergoing medical treatment, e.g. intra-vascular access or intravascular or therapy, or an arterio-venous (AV) shunt.
- Do not use the cuff on people who have undergone a mastectomy.
- Do not place the cuff over wounds as this may cause further injury.
- Place the cuff on your upper arm only. Do not place the cuff on other parts of the body.
- The blood pressure monitor can only be operated with batteries.
- To conserve the batteries, the blood pressure monitor switches off automatically if you do not press any buttons for 30 seconds.
- The device is only intended for the purpose described in these instructions for use. The manufacturer is not liable for damage resulting from improper or careless use.

Instructions for storage and maintenance

- The blood pressure monitor is made from precision and electronic components. The accuracy of the measured values and service life of the device depend on its careful handling:
 - Protect the device from impacts, humidity, dirt, marked temperature fluctuations and direct sunlight.
 - Do not drop the device.
 - Do not use the device in the vicinity of strong electromagnetic fields and keep it away from radio systems or mobile telephones.

- Only use the cuff included with the delivery or original replacement parts. Otherwise incorrect measured values will be recorded.
- We recommend that the batteries be removed if the device will not be used for a prolonged period of time.

Notes on handling batteries

- If your skin or eyes come into contact with battery fluid, rinse the affected areas with water and seek medical assistance.
-  **Choking hazard!** Small children may swallow and choke on batteries. Store the batteries out of the reach of small children.
- Observe the plus (+) and minus (-) polarity signs.
- If a battery has leaked, put on protective gloves and clean the battery compartment with a dry cloth.
- Protect batteries from excessive heat.
-  **Risk of explosion!** Never throw batteries into a fire.
- Do not charge or short-circuit batteries.
- If the device is not to be used for a relatively long period, take the batteries out of the battery compartment.
- Use identical or equivalent battery types only.
- Always replace all batteries at the same time.
- Do not use rechargeable batteries.
- Do not disassemble, split or crush the batteries.

Instructions for repairs and disposal

- Batteries do not belong in household waste. Please dispose of empty batteries at the collection points intended for this purpose.

- Do not open the device. Failure to comply will result in voiding of the warranty.
- Do not repair or adjust the device yourself. Proper operation can no longer be guaranteed in this case.
- Repairs must only be carried out by Customer Services or authorised suppliers. Before making a claim, please check the batteries first and replace them if necessary.
- For environmental reasons, do not dispose of the device in the household waste at the end of its useful life. Dispose of the unit at a suitable local collection or recycling point. Dispose of the device in accordance with EC Directive – WEEE (Waste Electrical and Electronic Equipment). If you have any questions, please contact the local authorities responsible for waste disposal.


Notes on electromagnetic compatibility

- The device is suitable for use in all environments listed in these instructions for use, including domestic environments.
- The use of the device may be limited in the presence of electromagnetic disturbances. This could result in issues such as error messages or the failure of the display/device.
- Avoid using this device directly next to other devices or stacked on top of other devices, as this could lead to faulty operation. If, however, it is necessary to use the device in the manner stated, this device as well as the other devices must be monitored to ensure they are working properly.
- The use of accessories other than those specified or provided by the manufacturer of this device can lead to an increase in electromagnetic emissions or a decrease in the device's electromagnetic immunity; this can result in faulty operation.


- Failure to comply with the above can impair the performance of the device.

3. Device description

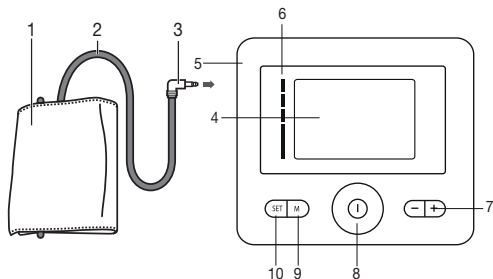1. Cuff
2. Cuff line
3. Cuff connector
4. Display
5. Connection for cuff connector
6. Risk indicator
7. Function buttons -/+
8. **START/STOP** button ①
9. Memory button **M**
10. Setting button **SET**

Information on the display:

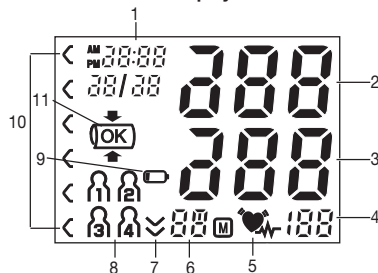
1. Time and date
2. Systolic pressure
3. Diastolic pressure
4. Calculated pulse value
5. Cardiac arrhythmia symbol
6. Memory display: average value (R), morning (RM), evening (PM), memory space number
7. Release air
8. User memory
9. Battery replacement symbol
10. Risk indicator
11. Cuff position control

4. Preparing the measurement

Inserting the batteries

- Open the battery compartment lid.
- Insert four 1.5V AA (alkaline type LR6) batteries. Make sure that the batteries are inserted the correct way round. Do not use rechargeable batteries.
- Close the battery compartment lid again carefully.
- All display elements are briefly displayed, 24 h flashes in the display. Now set the date and time as described below.

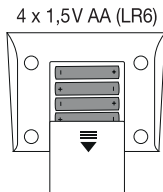If the low battery indicator  is permanently displayed, you can no longer perform any measurements and must replace the batteries. Once the batteries have been removed from the device, the time must be set again.

Battery disposal

- The empty, completely flat batteries must be disposed of through specially designated collection boxes, recycling points or electronics retailers. You are legally required to dispose of the batteries.
- The codes below are printed on batteries containing harmful substances:
Pb = Battery contains lead,
Cd = Battery contains cadmium,
Hg = Battery contains mercury

Setting the hour format, date and time

It is essential to set the date and time. Otherwise, you will not be able to save your measured values correctly with a date and time and access them again later.

There are two different ways to access the menu from which you can adjust the settings:

- Before initial use and after each time you replace the battery:
When inserting batteries into the device, you will be taken to the relevant menu automatically.
- If the batteries have already been inserted:
Press and hold the settings button **SET** on the device when **switched off** for approx. 5 seconds.

To set the date and time, proceed as follows:


- Insert the batteries, or press and hold the setting button **SET** for at least 5 seconds if the batteries have already been inserted.
- Select 12h or 24h mode using the function buttons -/+ . Press **SET** to confirm. The year display will start to flash. Set the year with the function buttons -/+ and confirm with **SET**.
- Set the month, day, hour and minute and confirm each with the setting button **SET**.
- The blood pressure monitor switches itself off automatically.

5. Measuring blood pressure


Ensure the device is at room temperature before measuring. The measurement can be performed on the left or right arm.

Attaching the cuff


Place the cuff onto the bare upper arm. The circulation of the arm must not be hindered by tight clothing or similar.

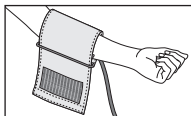
The cuff must be placed on the upper arm so that the bottom edge is positioned 2 – 3 cm above the elbow and over the artery. The line should point to the centre of the palm.

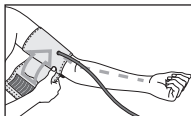
Now tighten the free end of the cuff, but make sure that it is not too tight around the arm and close the hook-and-loop fastener. The cuff should be fastened so that two fingers can fit under the cuff.

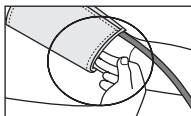Now insert the cuff line into the connection for the cuff connector.

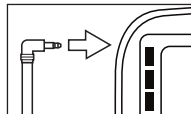This cuff is suitable for you if the index mark (▼) is within the OK range after fitting the cuff on the upper arm.



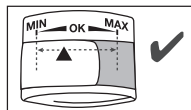
- i** If the measurement is performed on the right upper arm, the line should be located on the inside of your elbow. Ensure that your arm is not pressing on the line.

Blood pressure may vary between the right and left arm, which may mean that the measured blood pressure values are different. Always perform the measurement on the same arm. If the values between the two arms are significantly different, please consult your doctor to determine which arm should be used for the measurement.

Important: The unit may only be operated with the original cuff. The cuff is suitable for an arm circumference of 22 to 36 cm.

Adopting the correct posture

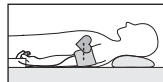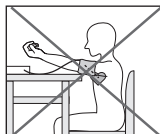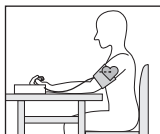- Before the initial blood pressure measurement, make sure always to rest for about 5 minutes. Otherwise deviations can occur.
- Furthermore, if you want to take several measurements in succession, make sure always to wait for at least 1 minute between the individual measurements.
- You can take the measurement while sitting or lying. Always make sure that the cuff is at heart level.
- To carry out a blood pressure measurement, make sure you are sitting comfortably with your arms and back leaning on

something. Do not cross your legs. Place your feet flat on the ground.

- To avoid falsifying the measurement, it is important to remain still during the measurement and not to speak.

Performing the blood pressure measurement

- As described above, attach the cuff and adopt the posture in which you want to perform the measurement.
- Start the blood pressure monitor with the **START/STOP** button ①.
- The cuff automatically inflates. The cuff's air pressure is slowly released. If you already recognise a tendency for high blood pressure, you should reinflate the cuff and increase the cuff's pressure again. As soon as a pulse is found, the pulse symbol ♥ will be displayed.
- The cuff position control symbol (OK) is displayed throughout the entire measurement. If the cuff is applied too tightly or too loosely, then (E) and "Er 3" are displayed. In such cases, the measurement is cancelled after approx. 5 seconds and the device switches itself off. Apply the cuff correctly and take a new measurement.
- The systolic pressure, diastolic pressure and pulse rate measurements are displayed.
- You can cancel the measurement at any time by pressing the **START/STOP** button ①.
- E_ appears if the measurement has not been performed properly. Observe the chapter on error messages/troubleshooting in these instructions for use and repeat the measurement.
- Now select the desired user memory by pressing the memory button M. If you do not select a user memory, the meas-

urement is stored in the most recently used user memory. The relevant 1, 2, 3 or 4 symbol appears on the display.

- To switch off, press the **START/STOP** button ①. If you forget to turn off the device, it will switch itself off automatically after approx. 3 minutes.

6. Evaluating results

Cardiac arrhythmia:

This device can identify potential disruptions of the heart rhythm when measuring and if necessary, indicates this after the measurement with the symbol ♥. This can be an indicator for arrhythmia. Arrhythmia is an illness in which the heart rhythm is abnormal because of flaws in the bioelectrical system that regulates the heartbeat. The symptoms (skipped or premature heart beats, pulse being slow or too fast) can be caused by factors such as heart disease, age, physical make-up, excess stimulants, stress or lack of sleep. Arrhythmia can only be determined through an examination by your doctor. If the symbol ♥ is shown on the display after the measurement has been taken, repeat the measurement. Please ensure that you rest for 5 minutes beforehand and do not speak or move during the measurement. If the symbol ♥ appears frequently, please consult your doctor. Self-diagnosis and treatment based on the measurements can be dangerous. Always follow your GP's instructions.

Risk indicator:

The measurements can be classified and evaluated in accordance with the following table.

However, these standard values serve only as a general guideline, as the individual blood pressure varies in different people and different age groups etc.

It is important to consult your doctor regularly for advice. Your doctor will tell you your individual values for normal blood pressure as well as the value above which your blood pressure is classified as dangerous.

The bar chart on the display and the scale on the unit show which category the recorded blood pressure values fall into. If the values of systole and diastole fall into two different categories (e.g. systole in the High normal category and diastole in the Normal category), the graphical classification on the device always shows the higher category; for the example given this would be High normal.

Blood pressure value category	Systole (in mmHg)	Diastole (in mmHg)	Action
Setting 3: severe hypertension	≥ 180	≥ 110	seek medical attention
Setting 2: moderate hypertension	160–179	100–109	seek medical attention
Setting 1: mild hypertension	140–159	90–99	regular monitoring by doctor
High normal	130–139	85–89	regular monitoring by doctor
Normal	120–129	80–84	self-monitoring

Blood pressure value category	Systole (in mmHg)	Diastole (in mmHg)	Action
Optimal	< 120	< 80	self-monitoring

Source: WHO, 1999 (World Health Organization)

Low blood pressure

Warning

Low blood pressure (hypotension) can be a health hazard and cause dizziness or fainting. Blood pressure is considered low if systolic and diastolic pressure are below 90/60 mmHG (source: National Health Service, 2023).

Seek medical attention if you suddenly suffer from low blood pressure.

7. Saving, displaying and deleting measured values

The results of every successful measurement are stored together with the date and time. If there are more than 30 measurements, the oldest measurements are lost.

- Press the memory button **M**. Select the desired user memory (M1 ... M4) by pressing the memory button **M** again.
- If you press the function button **+**, the average value $\bar{P}$ of all the stored measured values in the user memory will be displayed. If you press the function button **+** again, the average value of the morning measurements for the last 7 days will be displayed (morning: 5 a.m. – 9 a.m., display $\bar{P}_M^{\text{M}}$). If you press the function button **+** again, the average value of the evening measurements for the last 7 days will be displayed (evening: 6 p.m. – 8 p.m., display $\bar{P}_E^{\text{M}}$). If you continue to

press the function button **+**, the most recent individual measured values are displayed in turn with the date and time.

- To switch off, press the **START/STOP** button **⏻**.
- If you forget to switch off the device, it will switch itself off automatically after 30 seconds.
- If you want to delete the whole memory for a specific user, press the memory button **M**. Press and hold down the memory button **M** and the setting button **SET** simultaneously for 5 seconds.

8. Error messages/troubleshooting

In the event of errors, the error message **E_** appears on the screen.

Error messages may appear if:

- It was not possible to correctly record the pulse: **E 1**;
- No measurement could be taken: **E 2**;
- The cuff is fastened too tightly or loosely: **E 3**;
- Errors occur during the measurement: **E 4**;
- The pump pressure is higher than 300 mmHg: **E 5**;
- There is a system error. If this error message appears, please contact Customer Services: **E 6**.
- The batteries are almost empty: .

In such cases, repeat the measurement. Ensure that the cuff line is correctly attached and that you do not move or speak during the measurement. If necessary, reinsert or replace the batteries.

9. Cleaning and storing the device and cuff

- Clean the device and cuff carefully using a slightly damp cloth only.

- Do not use any cleaning agents or solvents.
- Under no circumstances hold the device and cuff under water, as this can cause liquid to enter and damage the device and cuff.
- If you store the device and cuff, do not place heavy objects on the device and cuff. Remove the batteries. The cuff line should not be bent sharply.

10. Technical specifications

Model No.	SBM 18
Measurement method	Oscillometric, non-invasive blood pressure measurement on the upper arm
Measurement range	Cuff pressure 0-300 mmHg, Systolic 50-280 mmHg, Diastolic 30-200 mmHg, Pulse 40-199 beats/minute
Display accuracy	Systolic ± 3 mmHg, Diastolic ± 3 mmHg, Pulse $\pm 5\%$ of the value shown
Measurement inaccuracy	Max. permissible standard deviation according to clinical testing: Systolic 8 mmHg / Diastolic 8 mmHg
Memory	4 x 30 memory spaces
Dimensions	L 110 mm x W 112 mm x H 53 mm
Weight	Approximately 341 g (without batteries, with cuff)
Cuff size	22 to 36 cm

Permissible operating conditions	+10°C to +40°C, 10-85% relative humidity (non-condensing)
Permissible storage conditions	-20°C to +55°C, 10-90% relative humidity, 800-1050 hPa ambient pressure
Power supply	4 x 1.5V — AA batteries
Battery life	For approx. 300 measurements, depending on levels of blood pressure and pump pressure
Classification	Internal supply, IP21, no AP or APG, continuous operation, application part type BF

The serial number is located on the device or in the battery compartment.

Technical information is subject to change without notification to allow for updates.

- This device complies with European Standard EN60601-1-2 (In accordance with CISPR 11, IEC 61000-4-2, IEC 61000-4-3 and IEC 61000-4-8) and is subject to particular precautions with regard to electromagnetic compatibility. Please note that portable and mobile HF communication systems may interfere with this unit. More details can be requested from the stated Customer Services address or found at the end of the instructions for use.
- The device complies with the EU Medical Devices Directive 93/42/EEC, the German Medical Devices Act (Medizinproduktegesetz) and the standards EN1060-3 (Non-invasive sphygmomanometers – Part 3: Supplementary requirements for electro-mechanical blood pressure measuring systems) and IEC80601-2-30 (Medical electrical equipment – Part 2 –

30: Particular requirements for the basic safety and essential performance of automated non-invasive sphygmomanometers).

- The accuracy of this blood pressure monitor has been carefully checked and developed with regard to a long useful life. If using the device for commercial medical purposes, it must be regularly tested for accuracy by appropriate means. Precise instructions for checking accuracy may be requested from the service address.

11. Replacement parts and wearing parts

Replacement parts and wearing parts are available from the corresponding listed service address under the stated material number.

Designation	Item number and/or order number
Standard cuff (22–36 cm)	163.912

12. Warranty / Service

Further information on the guarantee and guarantee conditions can be found in the guarantee leaflet supplied.

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